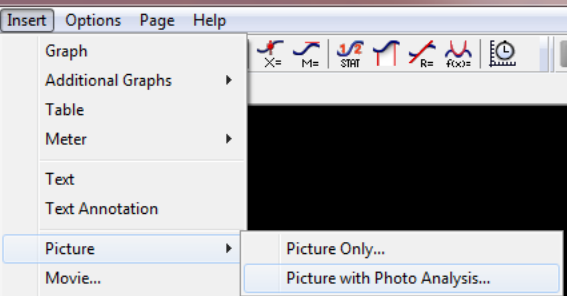
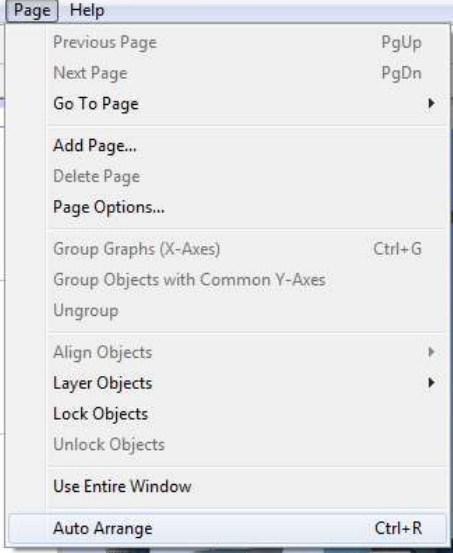







How to Use the Image Analysis Feature in Logger Pro

<ol style="list-style-type: none"> 1. Insert – Picture with Photo Analysis 2. Find image saved on your computer that you wish to use. 3. Resize picture by dragging corners. 	
<ol style="list-style-type: none"> 4. Auto Arrange Windows – this allows you to view graph, picture and table 5. Click on Page>Auto Arrange 	
<ol style="list-style-type: none"> 6. Find a measurement that you know is accurate. Ex. Width of St. Louis Arch is 630 ft. 7. Use Set Scale Tool  to draw a line on the measurement that you know. 8. Enter the measurement and unit in the box and Click OK. 	
<ol style="list-style-type: none"> 9. Then use the Set Origin Tool  to set the origin of your image. 10. Click on the picture where you want the origin after you have selected the tool. You will see the axes appear in yellow. 	

How to Use the Image Analysis Feature in Logger Pro

11. Use the Add Point Tool  to add data points to your image.
12. Click the image where you want the data points to appear. The data points will appear as red dots.

